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INSTRUCTIONS ON ENVIRONMENT PROTECTION

- At the end of its life cycle, please, do not dispose of this equipment by throwing it in the usual garbage; hand it over a recycling collection point.
- It does not contain dangerous or toxic products for humans but a non-adequate disposal would damage the environment.

- Materials are recyclable as mention in its marking.

- By recycling material or by other form of re-utilization of old appliances, you are making and important contribution to protect our environment.

- Please contact the administration of your community for being advised on the collection points.

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This manual should be available for all users of these equipments. To get the best results and a higher duration of this equipment it is advisable to read carefully this manual and follow the processes of use.

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f) Turn off the vacuum pump once sample filtration has finished and then disconnect the hoses and remove the filtration funnels. Please refer to “*Maintenance and cleaning*” point for cleaning all parts.

Security

- ♦ The filtration manifold must be used by previously qualified staff that knows how the equipment works thanks to the user manual.
- ♦ Establish all the security measures required for the occupational health and safety regulations in force (electronic devices, flammable, toxic or pathological products, etc.).
- ♦ Put the filtration manifold on top of a horizontal, plane and stable table making the necessary free space around.
- ♦ Do not put the filtration manifold near any warm supply (burners, blowlamps...), nor expose it directly to the sun, etc.
- ♦ Never dismantle the different pieces of the Nahita filtration manifold to repair it by yourself, since it could produce a defective use of the whole equipment and a loss of the product warranty, as well as injuries on people that handle the equipment.

6. MAINTENANCE AND CLEANING

All the processes of use mentioned below will not have any value unless you keep a continued and careful maintenance. To get the best results and a higher duration of this equipment it is essential to follow the processes of use.

- ♦ Please follow the processes of use of this manual.
- ♦ This manual should be available for all users of this equipment.
- ♦ Always use original components and supplies. Other devices can be similar but they can damage the equipment.
- ♦ Never try to repair the filtration manifold by yourself, since you will lose the warranty and may provoke damages to the general operating system of the equipment.
- ♦ In the event of breakdown please contact your distributor to overhaul through Nahita Technical Assistance Department.



Cleaning

- ♦ All components of manifolds and filtration holders can be sterilized in autoclave (121°C) or by dry heat.
- ♦ Never use scourers or substances that can scratch metallic parts such as stainless steel or aluminium, as they can damage the filter holder manifold and produce an early ageing of the equipment.
- ♦ Always use a sponge with warm soaped water for cleaning stainless steel or aluminium parts. You can also use products that do not contain abrasive or alkaline substances.


ATTENTION!! IF EQUIPMENTS ARE NOT PROPERLY CLEAN AND DISINFECTED THEY WOULD NOT BE ALLOWED TO REPAIR BY OUR TECHNICAL SERVICE.

b) Connect the manifold (*fig. 1.4*) and the filtering flask with a rubber hose and through the hole of the rubber stopper placed on the flask.

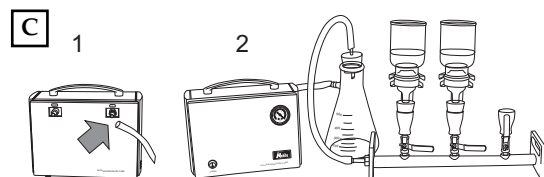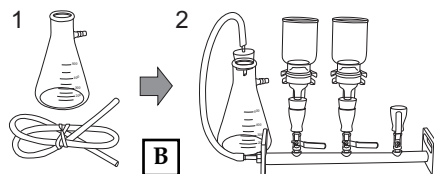c) Connect the filtering flask with the vacuum pump.

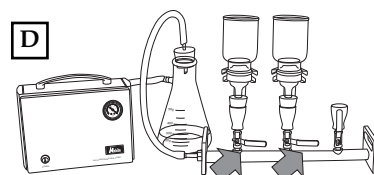
d) Pour the sample in the filtration funnel and make sure all the control valves (*Fig. 1.2*) you are working with are open.

e) Switch the vacuum pump on to begin sample filtration.

f) Turn off the vacuum pump once sample filtration has finished and then disconnect the hoses and remove the filtration funnel. Please refer to “Maintenance and cleaning” point for cleaning all parts.

Setting up, double circuit

a) Assemble the glass or stainless steel filter holder in each position according to point a) of Single Circuit manifolds. In case of using glass filter holders, don't forget to put a silicon lid on top of each funnel.

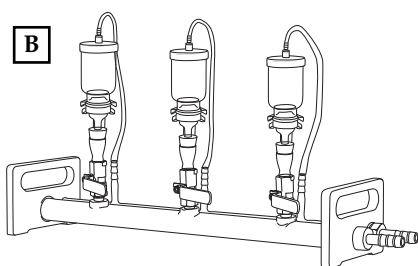b) Connect the hole at the top of each lid with the corresponding individual feeding branch so as the sample can reach funnel continuously (*fig. 2.5*).

c) Connect the tank with the liquid to be filtered with the feeding tube through the feeding connection (*fig. 2.7*). As an option and in case it is necessary a peristaltic pump can be used for supplying the sample (not included).

d) Connect the vacuum connection (*fig. 2.4*) of collector tube to a filtering flask and connect this flask to a vacuum pump, as previously shown in points b) and c) of Single Circuit manifolds.

e) Make sure all the control valves (*fig. 2.2*) you are working with are open and turn on the vacuum pump; this way the circuit of continuous feeding and filtration will be established. Note: always check the correct flow of sample.

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

Setting up, single circuit

a) Assemble the filter holder on each position of the manifold both for glass filter holder (fig.3) or for stainless steel filter holder (fig.4).

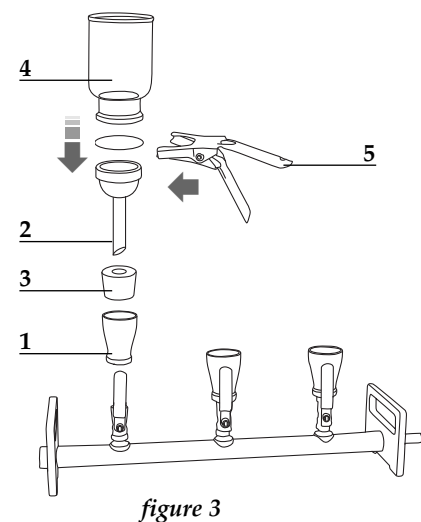- ◆ Screw the glass filter holder-manifold connector (1).
- ◆ Fasten properly the body with fritted (2) disk together with the silicone stopper (3) on the connector.
- ◆ Place the membrane filter.
- ◆ Place the filtration funnel (4) just on top of it and fasten the set with the aluminum clamp (5).

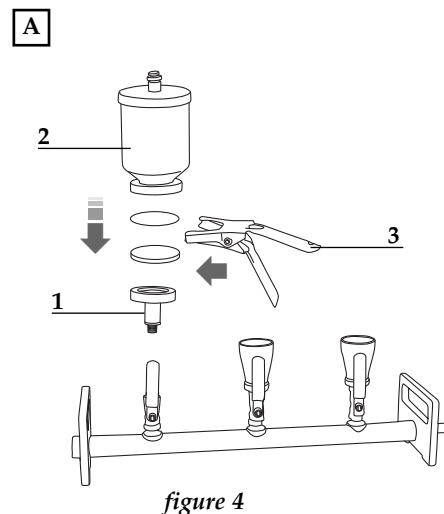- ◆ Screw the filtration body (1)
- ◆ Put the grating in its place on the filtration body
- ◆ Place the membrane filter
- ◆ Place the filtration funnel (2) just on top of it and fasten the set with the aluminium clamp (3).

NOTE: When screwing the connector (glass filter holder) or filtration body (stainless steel filter holder) to the manifold it is highly recommended to use PTFE sealing tape to assure a good sealing when making vacuum.

5. INSTALLATION / SETTING UP

Preliminary inspection

- ◆ Unwrap the filtration manifold, take off the involving plastic and take off the polyspan protection in which it comes fitted.
- ◆ Make sure that it does not present any damage because of the shipment. If so, tell it immediately to your transport agent or dealer so that they can make the claims in the correct time limit.

Please keep the original wrapping; you will always need it for returns enclosed with all the accessories supplied.

- ◆ Please check that all the accessories are enclosed with the equipment:
 - Filter holder manifold
 - Connection tubes
 - User manual

Installation

- ◆ Before using this instrument, it is convenient for you to familiarize with its components and basic essentials.

PLEASE READ THOROUGHLY THE INSTRUCTIONS BEFORE CONNECTING AND OPERATING WITH THIS EQUIPMENT

- ◆ Put the filtration manifold on top of a horizontal, plane and stable table making the necessary free space around.
- ◆ Establish all the security measures required for the occupational health and safety regulations in force (electronic devices, flammable, toxic or pathological products, etc.).
- ◆ Do not put the filtration manifold near any warm supply (burners, blowlamps...), nor expose it directly to the sun, etc.

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
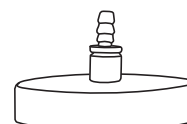
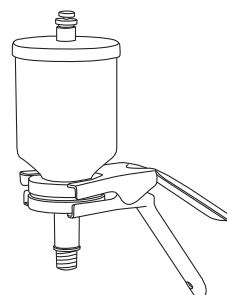
For a more detailed information please visit our website www.auxilab.com

4. ACCESSORIES (Not included)

The following glass and SS filter holders and accessories will complete manifold assembly according to user necessities.

Both filter holder sets available can be sterilized by autoclave (121°C) and dry heat.

Filter holders sets

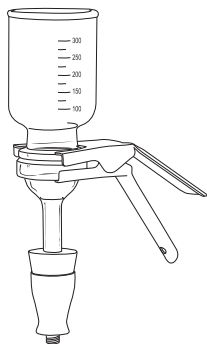Code	Description
15002900	Glass filter folder for manifolds


The following codes are included:

18002301	Filtration funnel, 300 mL
18003311	Body with fritted glass, 40 mm Ø
47195000	Silicon cap
65002300	Aluminium clamp
99050300	Glass filter holder-manifold connector

For using with double circuits manifolds it is also necessary to include the corresponding lids for glass funnels.

Code	Description
99051000	Silicon lid for double circuit manifold

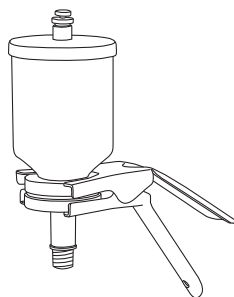Referencia	Descripción
99050330	Stainless steel filter holder 300ml
99050350	Stainless steel filter holder 500ml

The following codes are included:

99050303	Stainless steel filtration funnel 300 mL lid for double circuit manifold included (only in 99050330)
99050305	Stainless steel filtration funnel 500 mL lid for double circuit manifold included (only in 99050350)
99050310	Filtration body, stainless steel
99050311	Grating, stainless steel
65002300	Aluminium clamp

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
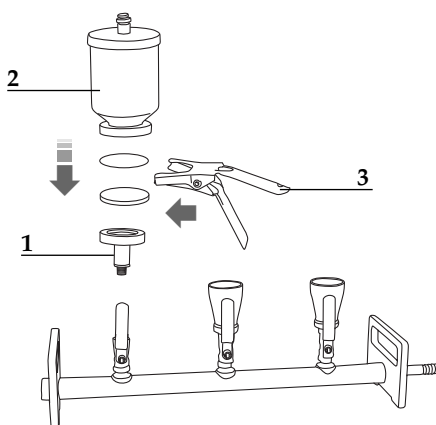
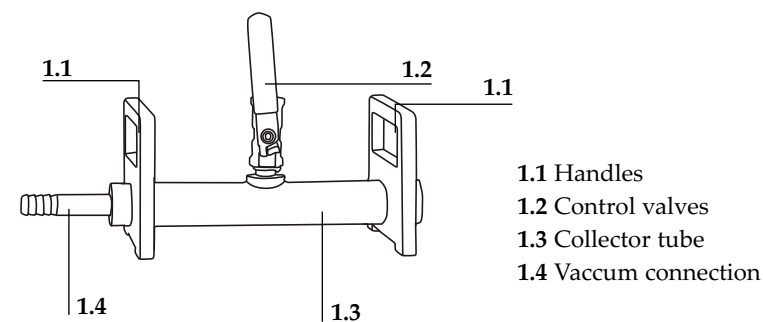
2. DESCRIPTION**Figure 1. Single circuit**

- 1.1 Handles
- 1.2 Control valves
- 1.3 Collector tube
- 1.4 Vacuum connection

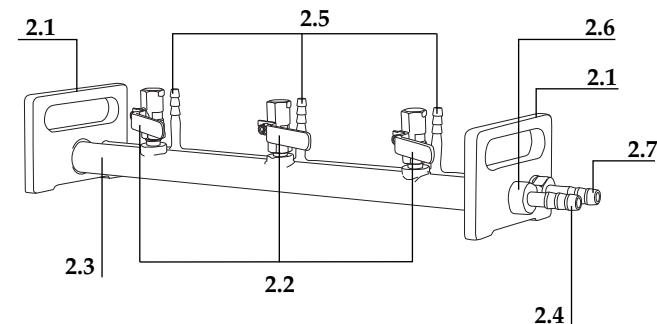Figure 2. Double circuit

- 2.1 Handles
- 2.2 Control valves
- 2.3 Collector tube
- 2.4 Vacuum connection
- 2.5 Individual feeding ranches
- 2.6 Feeding tube
- 2.7 Feeding connection

3. TECHNICAL SPECIFICATIONS

Code	Circuit	Positions	Body	Hose connection	Dimensions (mm)
59050301	Single	1	SS	12 mm	240x120x90
59050303	Single	3	SS	12 mm	480x120x90
59050306	Single	6	SS	12 mm	840x120x90
59051003	Double	3	SS	12 mm	470x120x90

INDEX OF LANGUAGES

Spanish.....	2-10
English.....	11-19

INDEX OF CONTENTS

1. USES OF THE INSTRUMENT.....	12
2. DESCRIPTION.....	13
3. TECHNICAL SPECIFICATIONS.....	13
4. ACCESSORIES.....	14
5. INSTALLATION / SETTING UP.....	16
6. MAINTENANCE AND CLEANING.....	19

1. USES OF THE INSTRUMENT

By using glass or stainless steel filter holders (**not included**), filtration manifolds allow vacuum filtration of one or several samples simultaneously. They have been specially designed to work in quality control of food, chemical or pharmaceutical industry as well as in water analysis.

The filtration manifold with double circuit (3 positions) allows a continuous filtration of big volumes of liquid, thus being ideal for microbiological control, control of pharmaceutical products and cosmetics, sterilization and preparation of big volumes of liquids and solvents, etc.

All available models are made of stainless steel, thus allowing sterilization by autoclave (121°C) or dry heat and present control valves (90°) in each position; so filtration in series with a common drainage can be made easy and quickly.

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Thank you for choosing this equipment. We sincerely wish that you enjoy your Nahita filter holder manifold. We highly recommend looking after this equipment according to what is stated in this manual.

The correct using of the equipment and its good quality will permit you to enjoy this equipment for years.

TO GET THE BEST RESULTS AND A HIGHER DURATION OF THE EQUIPMENT IT IS ADVISABLE TO READ THOROUGHLY THIS MANUAL BEFORE OPERATING WITH THE EQUIPMENT.

Please bear in mind the following:

- ♦ This manual is inseparable from the Nahita filter holder manifold, so it should be available for all the users of this equipment.
- ♦ You should carefully handle the Nahita filter holder manifold avoiding sudden movements, knocks, free fall of heavy / sharp objects on it.
- ♦ Never dismantle the different pieces of the Nahita filter holder manifold to repair it yourself, since it could produce a defective use of the whole equipment and a loss of the product warranty, as well as injuries on people that handle the equipment.
- ♦ If you have any doubt about setting up, installation or functioning do not hesitate in contacting your wholesaler.
- ♦ This equipment is protected under the Warranties and consumer goods regulation (10/2003).
- ♦ Overhaul is not covered by the Nahita warranty.
- ♦ Operations made by non-qualified staff will automatically produce a loss of the equipment's warranty.
- ♦ Accessories are not covered by the product warranty. The warranty neither covers piece's deterioration due to the course of time.
- ♦ Please make sure you keep the invoice, either for having the right to claim or asking for warranty coverage. In case you have to send the equipment to Nahita Technical Assistance Department you should enclose the original invoice or a copy as guarantee.
- ♦ Manufacturer reserves the right to modify or improve the manual or equipment.

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